



VISKOA®

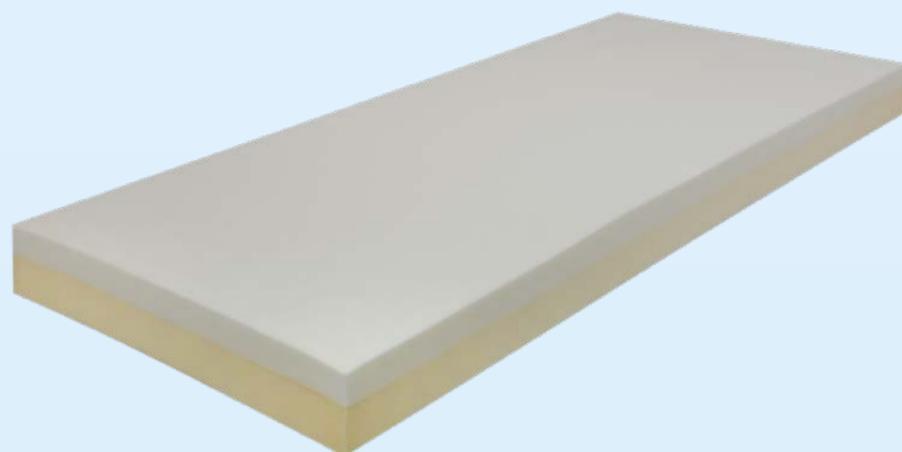

Braden scale score

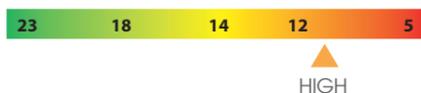

INDICATIONS

HELP IN PREVENTING PRESSURE ULCERS FOR PATIENT AT MEDIUM TO HIGH RISK OF PRESSURE ULCER DEVELOPMENT, BEDRIDDEN MORE THAN 15 HOURS A DAY.

HELP IN TREATING EXISTING PRESSURE ULCERS FROM STADE 1 TO 2 COMBINED WITH TECHNICAL AID DEVICES.

PATIENT WEIGHT: 30 TO 150 KG/ 4.7 TO 23.6 ST

- ▶ Single-piece mattress made up of visco elastic foam top layer on High Resilient (HR) polyurethane foam
- ▶ Fire Retardant properties in case of electric or fire ignition source
- ▶ Removable and washable integral cover



KEY FEATURES

- ▶ Single-piece mattress made up of:
 - 6 cm visco elastic top layer: density: 60 Kg/m³ ; hardness: 1.2 kPa
 - 8 cm High Resilient (HR) foam base: density: 35 Kg/m³ ; hardness: 3.9 kPa
- ▶ Patient weight: 30 to 150 Kg
- ▶ Dimensions: 195 x 88 x 14 cm (76.7 x 34.6 x 5.5 inches)
- ▶ Mattress weight: 14 Kg
- ▶ Cover available: Promust PU HD

KEY BENEFITS

- ▶ Visco elastic foam provides excellent pressure re-distribution properties: homogeneous immersion of the body into the foam, thus increasing the bearing surface and reducing pressure applied in areas at risk of pressure ulcers.
- ▶ Fire Resistant: complies with EN597-1 and 2, BS 6807: 2006 CRIB V and GPEM/CP D1-90: Classe B standards.
- ▶ Visco elastic foam provides patient comfort.

MATTRESS WITH COVER ⁽²⁾	TYPES OF IGNITION SOURCES THROUGH WHICH THE IGNITION RISK IS THE MOST IMPORTANT	PERFORMANCE CLASSES WITH RESPECT TO IGNITION ALLOWING TO REDUCE IDENTIFIED RISKS		
		SHORT STAY UNIT	LONG STAY UNIT	PSYCHIATRY UNIT
MEDICAL DEVICES ⁽¹⁾	Cigarette and match	D	C ⁽¹⁾	A ⁽¹⁾

(1) : There may be derogation from this performance when the person responsible for the establishment considers that the aspects of care, hygiene and comfort are preponderant on the fire risk.

(2) : Mattresses concerned by Decree No 95-292 and the Directive 93/42 on Medical Devices. Only the original cover allows the maintenance of safety requirements of Decree and CE marking.

ACCEPTANCE CRITERIA FOR FIRE TESTS



- ▶ Lower flame persistence 120 sec after burner removal.
- ▶ Limited deterioration of the mattress from the point of application of the burner.
- ▶ Combustion less one hour after the burner removal.



DESCRIPTION	ITEM REFERENCES ACCORDING TO PACKAGING		FIRE STANDARDS
	CARDBOX	PALLET BOX	
VISKOA® Promust PU HD cover	VMAFR-HIPH	VMAFR-HIPH BOX	EN 597-1:95 EN 597-2:95 GPEM/CP D1-90:1991

Cleaning instructions:
Promust PU HD cover



2 years warranty
1 year warranty for the cover in nursing homes, hospitals

Photos not contractually binding

