



ELK Lifting Cushion

The ELK is a lifting cushion designed to lift from the floor.

Easy to fit

- Lifting cushion is inflated using an Airflo Compressor with a simple hand control
- Simple to use, requiring very little training
- Suitable for use inside or out
- Can be used by one carer

Simple to maintain

- Clean using a proprietary non-abrasive mild liquid cleaner.
- Easy to deflate, fold away and store
- The material is designed to exceptionally high standards with extraordinary durability, strength & performance in the toughest of environments

Safe to use


- Lifts the frailest person up to 980 lbs from the floor
- Exceptionally rigid and stable
- Ideal for carers of individuals with neurological disorders who fall regularly

Improves well-being

- Provides a safe, dignified lift for a fallen person
- Reduces the risk of musculoskeletal injury for the carer
- Can be used in confined spaces where traditional lifting equipment may not fit

Technical Specifications

Max user weight	450 kg
Max height	56 cm
Min height	Flat
Seat depth	50 cm
Width	57"
Total product weight	3.6 kg
Airflo Compressor	24
Lifting performance per charge	Up to 9 lifts

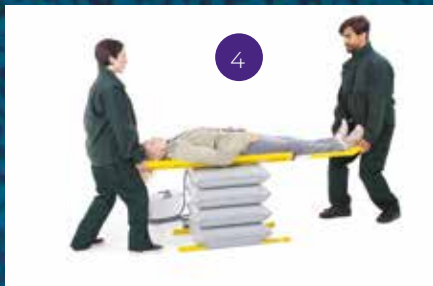

Product code: MPCA070400

Features

1. **4-Way Hand Control.**
2. **Made from very strong polyurethane** coated nylon material.
3. **Remove or keep stability bar** to provide extra comfort for the user.
4. **Can be used to lift a stretcher.**
5. **Easy to setup** - connect the four hoses from the hand control to the corresponding numbers and coloured connectors on the Camel.
6. **Powered by an Airflo 24 Compressor.**

Accessories

- **ELK Stretcher Bars**
- **Slidesheet**
- **ELK End Cap Set**



For more information visit

www.mangarhealth.com

Email: sales@mangarhealth.com

Phone: 0800 2800 485

