





Airflo Mk3		Airflo Plus	
<ul style="list-style-type: none"> • 12V - 'Normal' air flow only (No 'Normal - High' flow switch) • Battery only, no mains option • Plug top power supply (charger) - charging control is on the Airflo pcb • Integral LED battery indicator • No integral control buttons • 'Paddle' Hand Control included <p><i>Airflo Mk3 Supplied to operate with:</i></p> <ul style="list-style-type: none"> • Bathing Cushion • Handy Pillowlift • Sit-u-Up • Leglifter • Surfer Bather <p>Colour coded (Lilac) handle, label and power supply label for easy identification.</p>		<ul style="list-style-type: none"> • 24V - 'High' air flow only (No 'Normal - High' flow switch) • Battery and mains versions available • Plug top power supply (charger) - charging control is on the Airflo pcb • Integral LED battery indicator • Integral 'Auto' button with indicator LED • Hand Control NOT included <p><i>Airflo Plus Supplied to operate with:</i></p> <ul style="list-style-type: none"> • ELK • Camel <p>Colour coded (Light Grey) handle, label and power supply label for easy identification.</p>	
Battery: 12v 2.3Ah fixed (1 battery)		Battery: 24v 2.3Ah removable pack (2 batteries)	
Charging power supply: 15vDC 350mA with 2.1mm internal diameter DC output plug		Charging power supply: 30vDC 450mA with 2.5mm internal diameter DC output plug	
Unit weight: 4.6Kg	Pressure: 0.35bar (5psi)	Unit weight: 5.7Kg	Pressure: 0.35bar (5psi)
Noise: Noise level suitable for use in bedroom and bathroom.		Noise: Emits more noise than the Mk3 and so is less suitable for use in bedroom and bathroom.	
Lifting performance (approx from full battery charge):		Lifting performance (from full battery charge):	
Bathing Cushion	8 Lifts (4 Bathing cycles)	ELK	3 Lifts
Handy Pillow Lift	11 Lifts	Camel	1 Lift
Sit-u-Up	20 Lifts		
Leglifter	8 Lifts		
Surfer Bather	20 Lifts		
Charge Recovery time to do one further lift - 1Hr		Charge Recovery time to do one further lift - 1½ Hr	

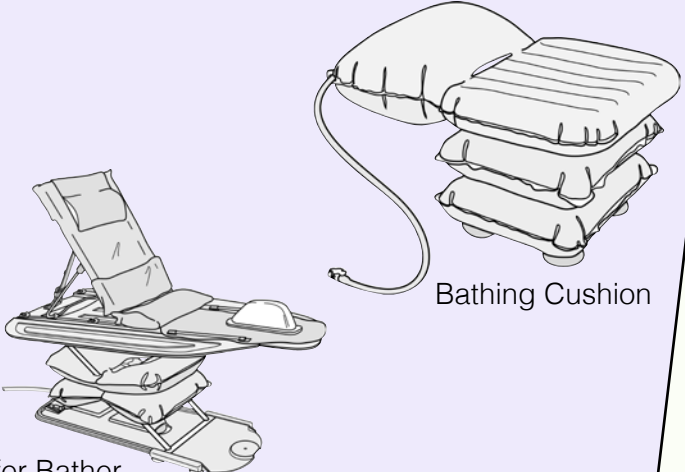

Surfer Bather


Bathing Cushion

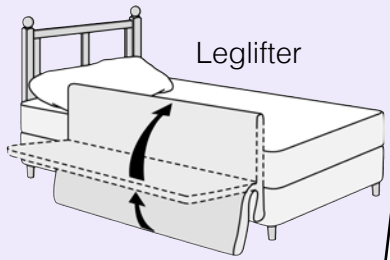


Leglifter




Sit-u-Up

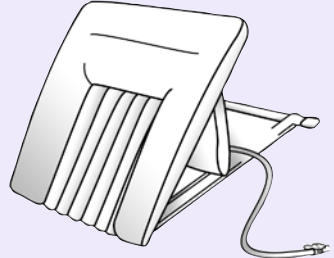


Handy Pillow Lift

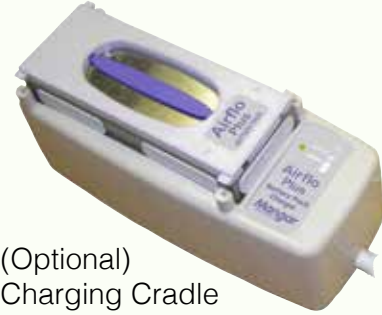

(Optional)
Charging Cradle

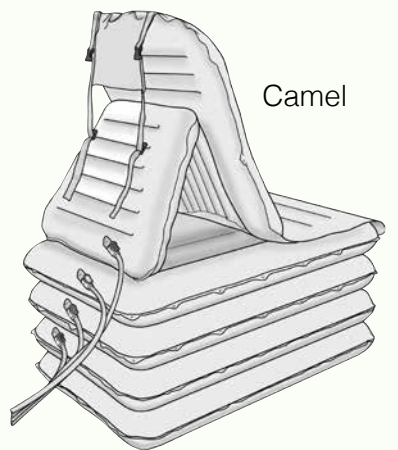

Camel

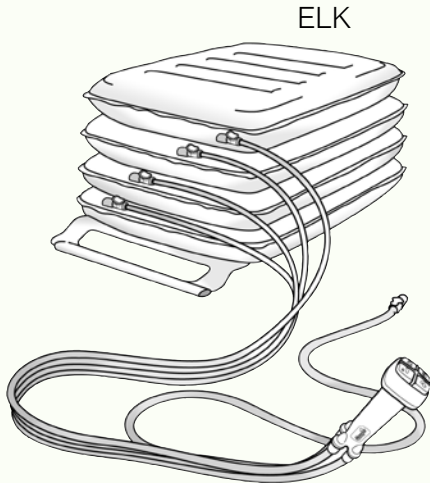

ELK